



No nonsense, all expertise

Our replaceable heads are designed by experts with patented bristle technology for a phenomenal clean.

GO THE DISTANCE

Help keep your mouth clean, replace your toothbrush head every 3 months.

YOU DO YOU

Find the perfect match for your hum brush in-store or online.

Apple®, App Store® and the Apple logo® are registered trademarks of Apple, Inc.

Google Play and the Google Play logo are trademarks of Google LLC.



Download on the
App Store

Saves information
to the Health app



GET IT ON
Google Play



powered by

Kolibree

DATA CHARGES MAY APPLY.

HERE TO HELP

1.833.hum.0001
(1.833.486.0001)

OR

Find hum online

Dist. by:
COLGATE-PALMOLIVE COMPANY
New York, NY 10022 U.S.A.
Made in China

Patent: www.colgatepalmolive.com/en-us/patents

P9941544 NAHUMB1SPLA1

READ ALL INSTRUCTIONS BEFORE USING

WARNING

For adults and children 8 years and older. Adult supervision required for persons with disabilities. Replace refill if brush head becomes loose due to damage. Battery choke hazard. Do not install battery backwards, charge, put in fire, or mix with other battery types. May explode or leak causing injury. Replace all batteries at the same time.

SAVE THESE INSTRUCTIONS

INSTRUCTIONS

CLEANING INSTRUCTIONS

After use, rinse toothbrush head thoroughly, shake off excess water, and store upright.

TO REPLACE BATTERIES

Firmly hold the bottom of the handle while rotating the top anti-clockwise. Insert 2 new AAA batteries as indicated inside the brush handle and firmly rotate clockwise back into position to close.

BLUETOOTH® 4.0

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and are used under license.

HANDLE: MODEL ADVTBC, BLUETOOTH® 4.0

BATTERY REMOVAL FOR RECYCLING

CORRECT DISPOSAL



Do not throw away this product with normal household waste at the end of its life. Always remove batteries before discarding. Please use a return collection system to ensure responsible and environmentally safe recycling of the device and battery.



REMOVE BATTERY BEFORE DISPOSAL

Remove batteries in compliance with environmental regulations before disposing of the toothbrush.

Ni-MH

FCC STATEMENT

Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help

The device must not be co-located or operating in conjunction with any other antenna or transmitter.

This device complies with Part 15 of the FCC Rules.

Operation is subject to the following two conditions :

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

WARRANTY

You can view the warranty of your purchased product by visiting our website: **www.colgate.com/humwarranty**

GET STARTED


1. INSERT BATTERIES

- Insert 2 new AAA batteries as indicated inside the brush handle



2. SET UP YOUR TOOTHBRUSH


- download the hum by Colgate app
- launch the app and check if Bluetooth® is turned on
- create your profile

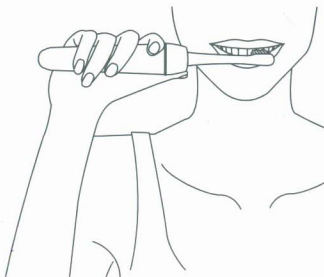

3. TURN YOUR TOOTHBRUSH ON

- a) Brush following the instructions from the app on your phone, OR
- b) Brush off line: *Timer alerts every 30 seconds to brush different areas of your mouth. Distinctive timer alert after 2 minutes. Turn the brush off when you're finished.*

Open the app and sync your toothbrush every 10 days to view your full brushing analysis and collect points toward rewards.

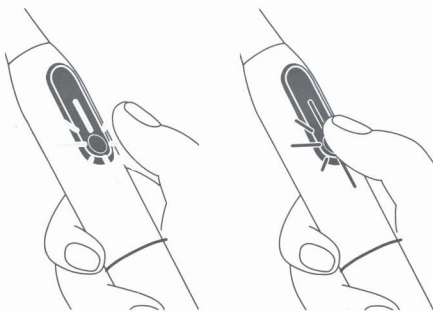

4. SWITCH BETWEEN MODES

Brush opens in the last cleaning mode used. Press button to switch to next mode: the LED indicates which mode you are in

- *White LED = NORMAL operates at high speed*
- *Pink LED = SENSITIVE operates at low speed*

Allow the brush to operate in a chosen mode for at least 3 seconds then press the power button to turn the brush off.